


Model DS8030 & DS8080

Wheelchair Scales



Model DS8030

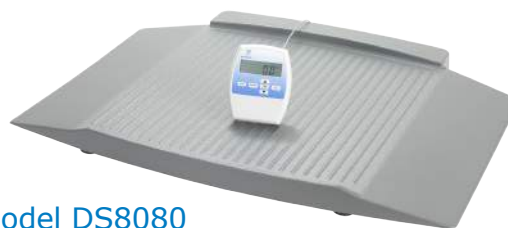

Model DS8080

The DS8030 and DS8080 offer a large, low profile base that makes weighing patients fast and easy. The weighing platform is constructed of heavy gauge painted steel and is sheathed with a textured, heavy-duty composite plastic cover that provides an anti-skid surface for maximum safety. The molded surface has no seams to trap dust and debris and is very easy to clean. The Model DS8080 features a digital readout that can be placed on a tabletop with its integrated stand, or mounted on a wall.

Two large diameter wheels are concealed along one side of these wheelchair scales to assist in moving the scales in your facility. With a few twists of the mast locking handles, the DS8030 folds down to a compact size so that it can easily be rolled through doorways. It can also be stored securely in this folded position thanks to the rubber stabilizer feet that are attached to the mast.

The DS8030 and DS8080 are EMR Ready and data can be transmitted via either the included USB or RS-232 ports.



MAXIMUM CAPACITY

800 lb x 0.2 lb (360 kg x 0.1 kg)

CONSTRUCTION

Painted Steel and Heavy Duty Composite Plastic

ACCURACY

+/- 0.2 lb (0.1 kg)

DIMENSIONS

DS8080	DS8030
Weighing Tray 30" W x 30" L	Weighing Tray 30" W x 30" L
Scale 37" W x 49" L x 3.5" H	Scale 37" W x 49" L x 50" H

DISPLAYED UNITS

lb or kg

POWER

AC Adapter or AA Batteries

WARRANTY

3 Three Years



CONTROLS

Units, Reweigh, Zero, Tare/BMI, Scroll Up/Down

