

Model DS9100


Wheelchair Scale

The ultimate in versatility and portability, the heavy-duty Model DS9100 Wheelchair Scale can serve as a wheelchair scale as well as a handrail or stand-on scale. To transport, the DS9100 folds and rolls effortlessly on four casters. Once in place, the DS9100 easily unfolds with built-in hydraulic assistance, and two latches hold the indicator mast firmly in place. The large, low-profile platform accommodates bariatric wheelchairs.

The DS9100 features a tare memory function that can store and recall the weights of up to ten different wheelchairs, eliminating error prone hand calculations to determine the patient's weight. To acquire a patient's weight, less the wheelchair tare weight, roll the patient onto the DS9100 platform, recall the wheelchair tare weight and an accurate weight up to 1,000 pounds will be held on the display. Our exclusive Weight Lock feature allows an accurate weight to be held regardless of a patient's motion.

The Reweigh feature allows a second weight to be acquired while the patient remains on the platform. Last Weight allows the previous weight and BMI value to be recalled after the patient has left the platform for later recording. To calculate BMI, enter the patient's height using the numeric keypad, press the BMI button and the patient's Body Mass Index is instantly displayed.

Optional USB and RS-232 communications options allow the recording of a patient's weight and BMI value into a PC based Excel spreadsheet with the optional Excelerator software. Patient data can also be sent to an optional printer or EMR system.

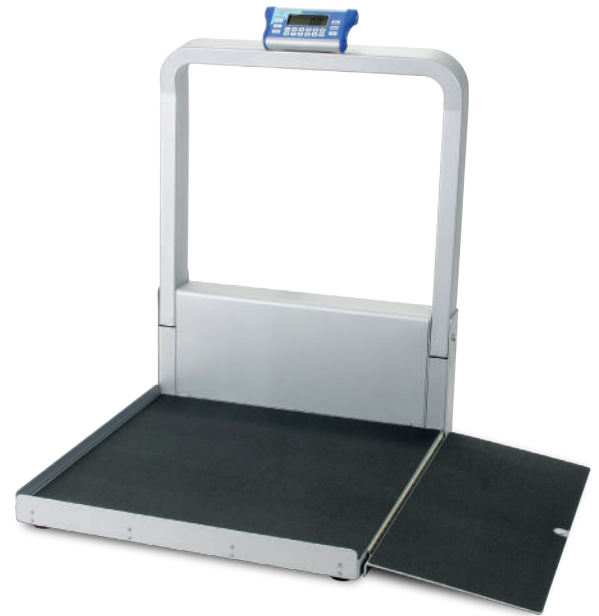

MAXIMUM CAPACITY

1000 x 0.1 lb (454 x 0.05 kg)

ACCURACY

+/- 0.1 lb (0.05 kg)

PLATFORM SIZE

33" W x 32.5" D x 3" H

SCALE DIMENSIONS

RAMP EXTENDED

49" W x 42" D x 45.5" H

DISPLAYED UNITS

lb, kg

DISPLAY

1" High LCD

BMI CALCULATOR

Input data in FT/IN or CM

APPROVALS

CE, UL (AC Adapter)

POWER

Six AA Batteries (Included)

BATTERY LIFE

One Year

CONSTRUCTION

Aluminum and Painted Steel
Construction with Plastic Trim

OPTIONS

(P/N DS9100-232)

RS-232 Communications
Port

(P/N DS9100-USB)

USB Communications Port
(P/N DS9100-RAMP)

Second Ramp

AC Adapter (CE or UL
Versions Available)

WARRANTY



Three Years



EMR Ready



Doran Scales, Inc. 1-800-264-4106

1315 Paramount Parkway Fax: 630-879-0073

Batavia, IL 60510 www.doranmedical.com