

ENOCHS HORIZON POWER 2250

High-Low Treatment Table

PATIENT LOAD: 500 LBS.
evenly distributed on the table's main couch

NAUGAHYDE® brand material with **ADVANCE BEAUTY GUARD®** protective finish is manufactured to have high resistance to abrasion, color fading, oil, mildew, sulfide stains and has agents effective against bacterial and fungal micro-organisms

ONE-PIECE UPHOLSTERED TOP
eliminates cross seams and decorative embossing to promote quick, efficient cleaning after each patient



The **PAPER ROLL HOLDER**
accommodates 18" or 21" wide paper

The **32" EXTRA WIDE TOP**
padded with 1.5" of high
density foam for patient's
comfort

'EASY TOUCH' FOOT CONTROL
adjusts table's height with reliable
hands-free action. Pedal is covered
in non-slip tread and utilizes the
natural left to right rocking motion
of the operator's foot.



The height adjustment range of 16" allows physician or staff member to position patient at a convenient working height. The Horizon 2250 is an excellent table for exam, procedural, treatment rooms as well as physical therapy and casting areas. The Easy Touch SELV Foot Control allows the operator to use the natural side-to-side movement of the foot to adjust the table's height from 23" to 39". Use of a foot control reduces incidents of cross contamination and enhances infection control procedures. Control pedal is covered with easy to clean non-slip tread.

The Horizon 2250 High-Low Treatment Table is an ideal choice for any medical facility requiring low access for their patients. The table's adjustable height and the 32" extra wide upholstered top are two distinctive features designed for patient comfort as well as physician and staff convenience. Table has a 500 pound patient load and offers easy low access by disabled or geriatric persons in wheelchairs or with walking aids, obese patients, and those requiring physical therapy. The low 23" access height helps reduce back stress or injury to staff assisting patients on and off table.

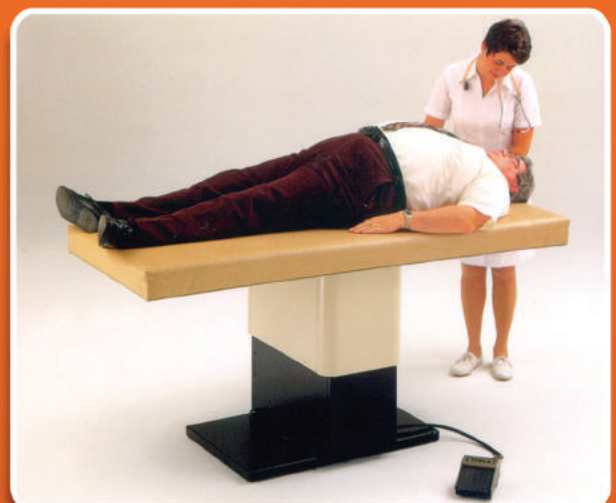


TABLE FEATURES:

- 4-YEAR Warranty
- Fully welded All Steel Base Construction - NO PLASTIC
- Full Power adjustment of height
- Easy to clean One-Piece Upholstered Top with 1.5" high density foam
- 32" Extra Wide Top for patient comfort
- SELV Foot Control
- Certified by UL to Standard 60601-1, CAN/CSA C22.3 No. 601-1
- ENOCHS Horizon 2250 replaces the Power 2000

Specifications for the ENOCHS Horizon 2250 High-Low Treatment Table

Warranty:	4-YEAR
Base Construction:	Fully welded all steel base
Power System:	Electro-Mechanical provides smooth and quiet movement
Power Adjustments:	Height: 23" - 39"
Foot Control:	"Hands free", standard on every table, SELV
Hand Control:	Not available on this model
Couch Top:	One-piece upholstered top
Perineal Recess:	Not available on this model
Pelvic Tilt:	Not available on this model
Leg Extension:	Not available on this model
Irrigation Basin:	Not available on this model
Electrical Outlets: with Circuit Breaker	Not available on this model
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model
Stirrups:	Not available on this model
Storage:	Not available on this model Total Storage Space: n/a
Paper Roll Holder:	Accommodates 18" or 21" paper widths
Accessories:	Not available for this model
Mains Power Entry: with Circuit Breaker	IEC connection, internal circuit breaker with reset
Duty Cycle:	Table: 1 minute ON / 9 minutes OFF Receptacles: n/a Instrument Warming Drawer: n/a
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: n/a Total Receptacle Load: n/a
Table Dimensions:	Width: 32", Height: 23" - 39", Length: 72"
Table Weight:	350 lbs.
Patient Load:	500 lbs. evenly distributed on main couch
Certifications:	By UL to Standard 60601-1, CAN/CSA C22.3 No. 601-1