

ENOCHS REGENCY 456

OB/GYN Examination Table

PATIENT LOAD: 600 LBS.
evenly distributed on the table's main couch

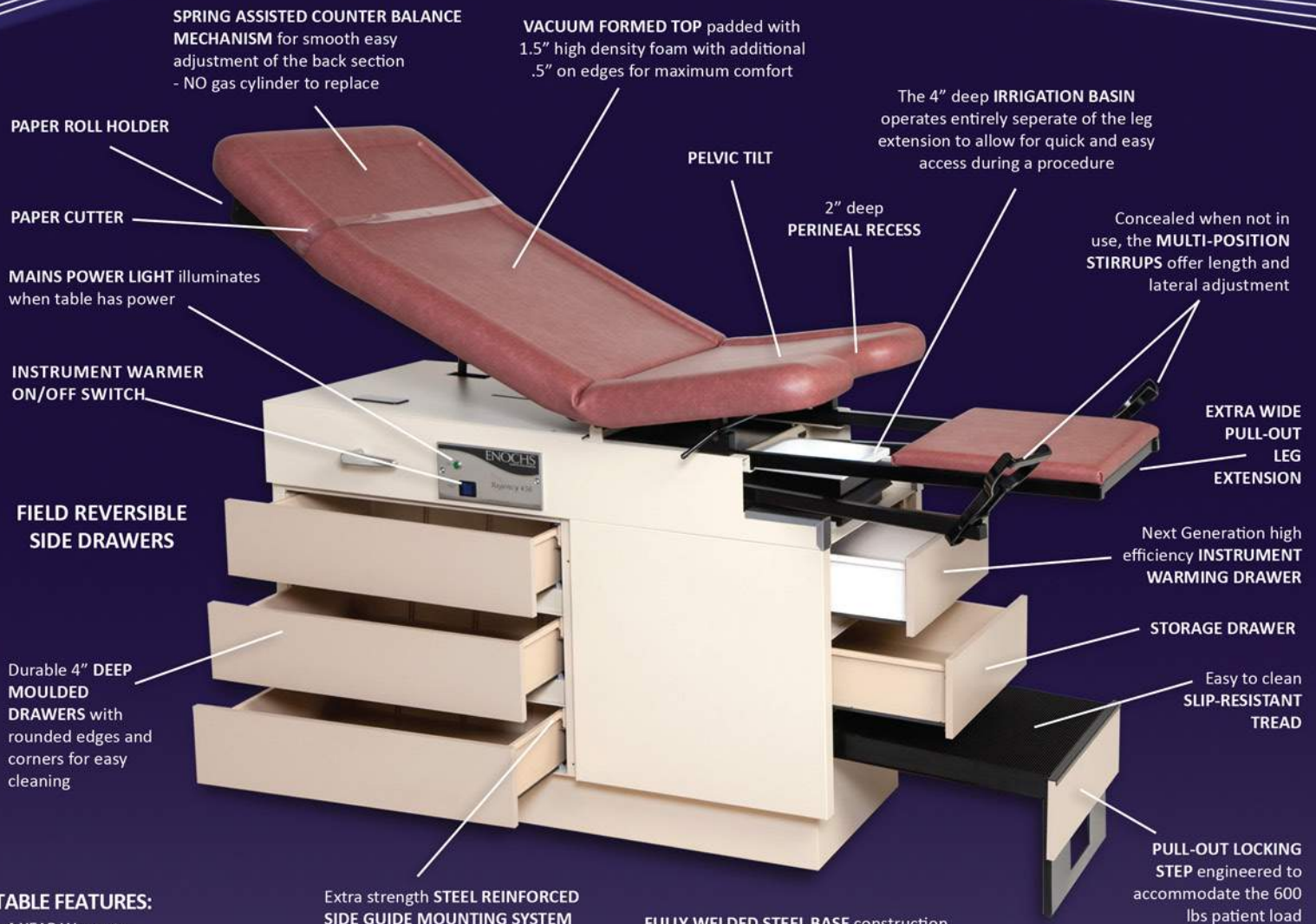


TABLE FEATURES:

- 4-YEAR Warranty
- Dual Electrical Outlets, Hospital Grade, 9 amps total load
- Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- ON/OFF Instrument Warmer Switch with integrated Circuit Breaker - NO FUSES
- All Steel Base and Top Construction - NO PLASTIC
- Easy to clean 100% Seamless Vacuum Formed Top with extra padding on edges for patient comfort
- Pelvic Tilt 7° & 11°
- Irrigation Basin designed to move independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers
- Certified by UL to Standard 60601-1 CAN/CSA C22.3 No. 601-1
- ENOCHS Regency 456 replaces the Regency 60

Extra strength **STEEL REINFORCED SIDE GUIDE MOUNTING SYSTEM**

FULLY WELDED STEEL BASE construction
- NO PLASTIC PARTS



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



Located at footend of table, the high efficiency instrument warmer is designed to warm speculums and instruments to approximate body temperature. The ON/OFF switch for the warmer is located on the front side of table for instant confirmation of unit's status. An internal circuit breaker in the switch will trip if circuit is overloaded, interrupting power to the warming unit but maintaining power to the electrical receptacles.

Specifications for the ENOCHS REGENCY 456 OB/GYN Examination Table

Warranty:	4-YEAR
Base Construction:	Fully welded all steel base
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance
Couch Top:	100% Seamless vacuum formed style All steel frame construction
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point
Pelvic Tilt:	7°, 11°
Leg Extension:	Extends out 20" from end of seat section
Irrigation Basin:	4" Deep basin
Electrical Outlets:	Dual electrical outlets, hospital grade, UL Listed
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	High efficiency warmer utilizes top footend drawer
Stirrups:	Multi-positional, lateral and length adjustment
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet
Foot Step:	Extra long surface engineered for 600 lb. patient load
Paper Roll Holder:	Accommodates 18" or 21" paper widths
Accessories:	knee crutches, arm board, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, headend accessory rails, footend accessory rails, caster base
Mains Power Entry: with Circuit Breaker	Standard feature
Duty Cycle:	Receptacles: Continuous Instrument Warming Drawer: Continuous
Ratings:	Table: 115V 60Hz 10 amps Receptacles: 115V 60Hz Total Receptacle Load: 9 amps
Table Dimensions:	Width: 27", Height: 34", Length: 52" (extends to 72")
Table Weight:	268 lbs.
Patient Load:	600 lbs. evenly distributed on main couch
Certifications:	By UL to Standard 60601-1, CAN/CSA C22.3 No. 601-1