

Enthermics Medical Systems


Warming patients inside
and out



 enthermics®

enthermics.com | ISO 13485:2003 certified company

Why Warm

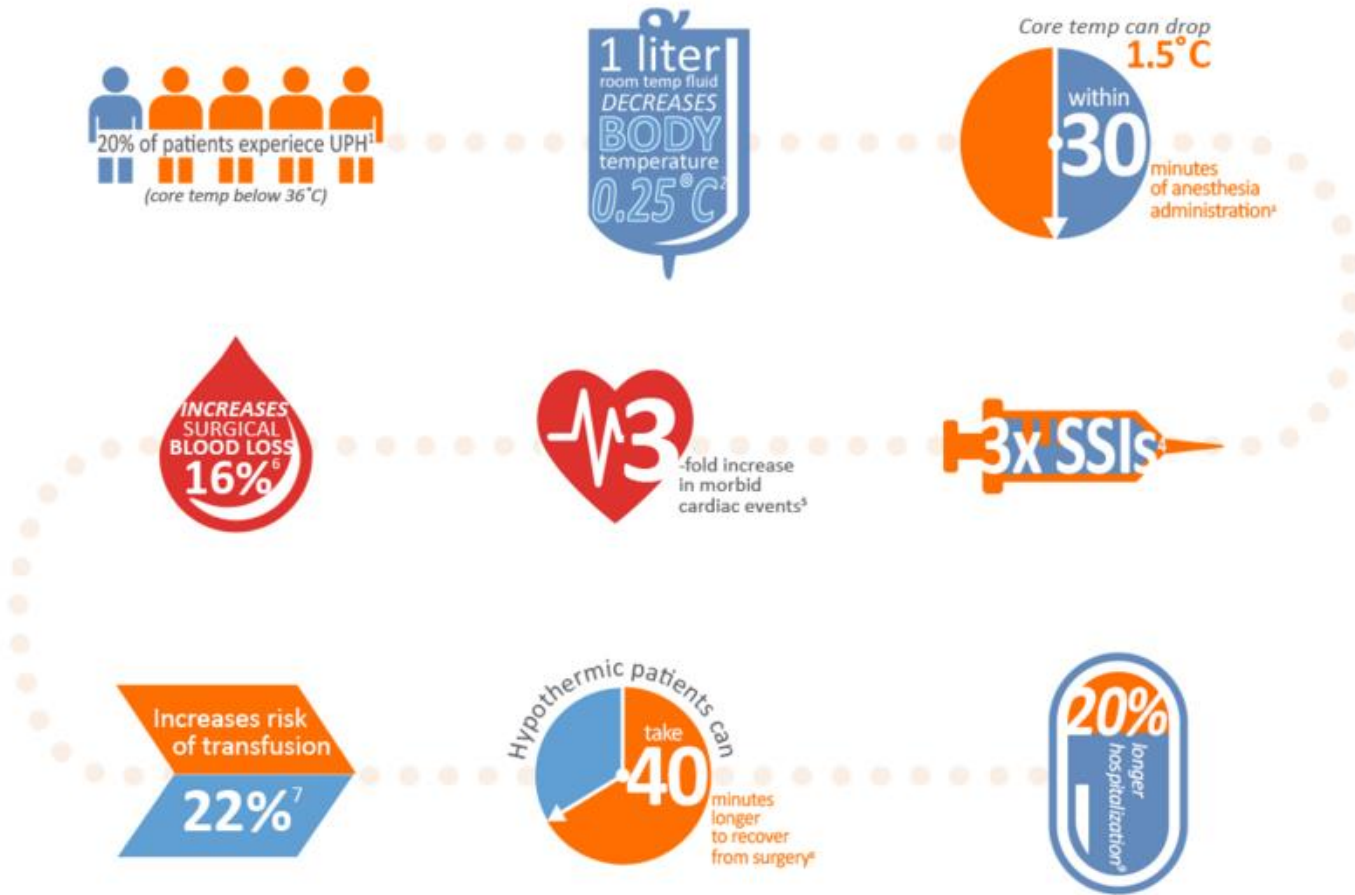

Methods of Patient Warming

Warm IV Fluids

Conductive Warming Devices

Warming Blankets

Forced Air

Circulating Water Garments

Increased Operating Room Temperature



Titan Series



- Energy efficient WarmSafe™ heating system monitors interior surface temperatures and balance heat throughout cabinet
(Patent No: US 8,217,316; US 8,581,152)
- Silent operation (no fans or moving parts)
- Lint traps not required
- Stationary or mobile option
- Epoxy-coated interior surfaces
- Programmable control
- Three capacities available:
 - » 4-6 blankets (2.5ft³)
 - » 10-12 blankets (3.5ft³)
 - » 20-25 blankets (7.5ft³)

Designer Series



- Energy-efficient WarmSafe™ heating system monitors interior surface temperatures and balances heat throughout the cabinet
(Patent No: US 8,217,316; US 8,581,152)
- Silent operation (no fans or moving parts)
- Lint traps not required
- Stationary or mobile option
- Epoxy-coated exterior & interior surfaces
- Programmable control
- Interior lighting with a comforting blue glow. Two (2) intensity settings
- Five capacities available:
 - » 3-4 blankets (1.5ft³)
 - » 4-6 blankets (2.5ft³)
 - » 10-12 blankets (3.5ft³)
 - » 14-18 blankets (4.0ft³)
 - » 20-25 blankets (7.5ft³)

No Fans in Blanket Warmers

- Hospitals do not get the chance to buy warmers very often.
- Sell - total cost of ownership
- No fans, no moving parts = Fewer repairs
- Broad range of sizes is a competitive advantage. Don't let it bog you down!

WarmSafe Technology

- Intelligent heating system prevents overheating of blankets
- System monitors temperature in each cabinet zone and provides heat when it is needed, where it is needed
- Overheating, blanket scorching and discoloration are prevented



Fluid Warmers

- Patented heating system
- Air mixing fan ensures accurate set-point temperature
- Fluids kept at proper temperature
- Eliminates disposable expense and set-up time of in-line warming systems



EC1850BL Combination Warmer



- Fluid Chamber holds 30 bottles of 1 liter or 24 bags of 1 liter.
- Blanket chamber capacity: 30 blankets (0.46 m³).
- WarmSafe™ heating system incorporates a multiple zone warming technology that heats where and when needed.
- Stainless steel exterior and interior.
- Glass doors
- Accurate temperature: IRR +0/-1.67°C; INJ +0/-1.1°C.
- Interior illumination
- 11.875 ft³ (0.34m³) chamber for blankets
- 38 liter fluid chamber capacity
- WarmRight® provides the right temperature for the right product: injection or irrigation fluids
- Blanket chamber: 32-71°C (90-160°F) range.
- Fluid chamber: Injection (INJ) or irrigation (IRR) fluid warming modes:
 - 32-40°C (90-104°F) adjustable INJ mode
 - 32-66°C (90-150°F) adjustable IRR mode

Quality

- ISO 13485:2003 certified company
- All fluid warmers are FDA cleared



Enthermics Quality Policy

Enthermics Medical Systems is committed to meeting our customers' expectations through continual improvement and providing high quality, innovative products with outstanding service from a company you can trust.



**All Enthermics products are
made in the USA**



enthermics.com | ISO 13485:2003 certified company